

## **Water Quality Monitoring: Total Dissolved Oxygen Assessment in Different Water Samples by Winkler's Method**

**1. Importance of Dissolved Oxygen (DO):** Knowledge of DO concentration in seawater is often necessary in environmental and marine sciences. It is used by oceanographers to study water masses in the ocean. It provides the marine biologist a means to measure primary production, particularly in laboratory cultures. For the marine chemist, it provides a measure of the redox potential of the water column. DO is also an important factor in corrosion. Oxygen is poorly soluble in water. The solubility of oxygen decreases with increase in concentration of the salt and hence, solubility of DO is lesser in saline water. The amount of DO at 100% saturation at sea level is 9.03 mg/L (at 20° C) and is sufficient to sustain aquatic life. Dissolved oxygen is usually determined by Winkler's method.

**2. What is Winkler Method?** The Winkler Method is a technique used to measure dissolved oxygen in freshwater systems. DO is used as an indicator of the water body's health, where higher DO concentrations are correlated with high productivity and little pollution. This test is performed on-site, as delays between sample collections and testing may result in an alteration in oxygen content.

**3. How does the Winkler Method Work?** Winkler Method uses titration to determine dissolved oxygen in the water sample. A sample bottle is filled completely with water (no air is left to skew the results). DO in the sample is then "fixed" by adding a series of reagents that form an acid compound that is then titrated with a neutralizing compound that results in a colour change. The point of colour change is called the "endpoint," which coincides with the dissolved oxygen concentration in the sample. DO analysis is best done in the field, as the sample will be less altered by atmospheric equilibration.

### **4. Applications:**


Dissolved oxygen analysis can be used to determine the health or cleanliness of a lake or stream, amount and type of biomass a freshwater, the amount of DO that a system can support and the amount of decomposition occurring in the lake or stream.

**Expt. No.:**

**Date:**

<b>Experiment</b>	<b>Water Quality Monitoring: Total Dissolved Oxygen Assessment in Different Water Samples by Winkler's Method</b>
<b>Problem definition</b>	Dissolved oxygen (DO) is essential to living organisms in water but harmful if present in boiler feed water leading to boiler corrosion.
<b>Methodology</b>	Winkler's titration method is used to assess DO in water.
<b>Solution</b>	Estimation of total dissolved oxygen in different water samples.
<b>Student learning outcomes</b>	Students will learn to a) perform Winkler's titration method b) assess the total dissolved oxygen in different water samples

**Principle:** Estimation of dissolved oxygen (DO) in water is useful in studying corrosion effect of boiler feed water and in studying water pollution. DO is usually determined by Winkler's titration method. It is based on the fact that DO oxidize potassium iodide (KI) to iodine. The liberated iodine is titrated against standard sodium thiosulphate solution using starch indicator. Since DO in water is in molecular state, as such it cannot oxidize KI. Hence, manganese hydroxide is used as an oxygen carrier to bring about the reaction between KI and Oxygen. Manganese hydroxide, in turn, is obtained by the action of NaOH on  $MnSO_4$ .

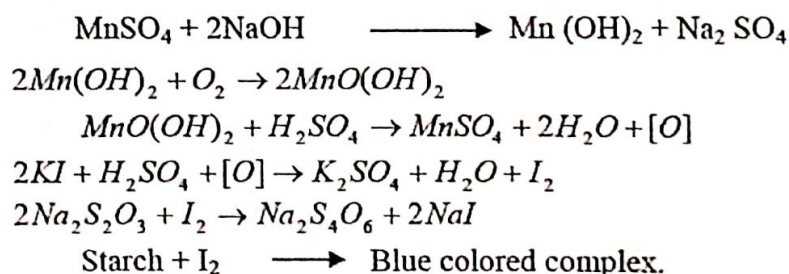

The liberated iodine ( $I_2$ ) is titrated against standard sodium thiosulphate ( $Na_2S_2O_3$ ) solution using starch as indicator.

**Requirements:**

**Reagents and solutions:** Standard buffer of pH 7, standard potassium dichromate (0.01 N), sodium thiosulphate solution, 10% KI solution, alkali KI solution (KI + KOH in water), conc.  $H_2SO_4$ , manganese sulphate, starch solution as indicator.

**Apparatus:** Conical flask, Burette, Measuring flask, Beakers.



### Procedure:

#### Titration 1: Standardization of Sodium Thiosulphate

Rinse and fill the burette with given sodium thiosulphate solution (Bottle B). Pipette out 20 mL of 0.01N  $K_2Cr_2O_7$  solution (Bottle A) into a clean conical flask. To this, add 5 mL  $H_2SO_4$  (1/2 T.T.), 5 mL (1/2 T.T.) of 10% KI, and titrate against sodium thiosulphate solution. When the solution becomes straw yellow colour, add starch indicator and continue the titration. End point is the disappearance of bluish brown colour. Repeat the titration to get concordant value.

#### Titration 2: Estimation of Dissolved Oxygen

Using a measuring cylinder, add 100 mL of sample water in a conical flask. Further, add 2 mL of  $MnSO_4$  and 2 mL of alkali KI solution and shake well for the rough mixing of the reagents. Set aside the flask for few minutes to allow the precipitate to settle down and then add 2 mL of conc.  $H_2SO_4$  for complete dissolution of the precipitate. Then, titrate against std. sodium thiosulphate solution. When the solution turn into light yellow, add starch indicator. End point is the disappearance of bluish brown colour. Repeat the titration to get the concordant value. Calculate the strength of dissolved oxygen from the titer value. Based on that, calculate the amount of dissolved oxygen in the given water sample.

### OBSERVATION AND CALCULATIONS

#### Titration - I: Standardization of Sodium Thiosulphate

S. No.	Volume of $K_2Cr_2O_7$ (mL)	Burette reading (mL)		Volume of sodium thiosulphate (mL)
		Initial	Final	
1	20	0	16.6	16.6
2	20	0	15.9	15.9
3	20	0	16.2	16.2
Concordant value				16.2

#### Calculations:

Volume of potassium dichromate  $V_1 = 20$  mL

Strength of potassium dichromate  $N_1 = 0.01$  N

Volume of sodium thiosulphate  $V_2 = 16.2$  mL (From Titration - I)

Strength of sodium thiosulphate  $N_2 = \dots\dots\dots?$

$$V_1 N_1 = V_2 N_2$$

$$\therefore N_2 = V_1 N_1 / V_2$$

Strength of sodium thiosulphate  $= N_2 = 20 \times 0.01 / 16.2 = 0.012$  N

**Titration – II: Estimation of Dissolved Oxygen**

S. No.	Volume of water sample (mL)	Burette reading (mL.)		Volume of sodium thiosulphate (mL)
		Initial	Final	
1.	100	0	6.5	6.5
2.	100	0	6.8	6.8
3.	100	0	6.5	6.5
Concordant value				6.6

**Calculation:**

Volume of sodium thiosulphate  $V_2 = 6.6$  mL (From Titration – 2)

Strength of sodium thiosulphate  $N_2 = 0.012$  N (From Titration – 1 calculation)

Volume of water sample taken  $V_1 = 100$  mL

Strength of given water sample  $N_1 = ?$

$$V_1 N_1 = V_2 N_2$$

$$N_1 = V_2 \times N_2 / 100$$

$$= 0.00072 \text{ N}$$

Amount of dissolved oxygen (ppm) = normality  $\times$  equivalent weight of  $O_2 \times 1000$  mg/L of the given water sample.

$$= 0.00072 \text{ N} \times 8 \times 1000 \text{ mg/L}$$

$$= 6.336 \text{ ppm.}$$

**Result:** Amount of dissolved oxygen in the given water sample = 6.336 ppm.

**Evaluation of Result:**

Sample number	Experimental value	Actual Value	Percentage of error	Marks awarded
			0.51. error	10/10